abcam

Product datasheet

Biotin Anti-AIF antibody [E20] - Mitochondrial Marker ab197526



重组 RabMAb

2 图像

概述

产品名称 生物素Anti-AlF抗体[E20] - Mitochondrial Marker

描述 生物素兔单克隆抗体[E20] to AIF - Mitochondrial Marker

宿主 Rabbit 偶联物 **Biotin**

经测试应用 适用于: IHC-P 种属反应性 与反应: Human

预测可用于: Mouse, Rat 📤

免疫原 Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.

阳性对照 IHC/P: Normal human colon tissue.

Our $\mathsf{RabMAb}^{\mathsf{®}}$ technology is a patented hybridoma-based technology for making rabbit 常规说明

monoclonal antibodies. For details on our patents, please refer to **RabMAb® patents**.

性能

形式 Liquid

存放说明 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C.

Avoid freeze / thaw cycle. Store In the Dark.

存储溶液 pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: 30% Glycerol (glycerin, glycerine), 1% BSA, PBS

纯度 Protein A purified

克隆 单克隆 克隆编号 E20 同种型 ΙgG

应用

Abpromise™承诺保证使用ab197526于以下的经测试应用 The Abpromise guarantee

"应用说明"部分 下显示的仅为推荐的起始稀释度:实际最佳的稀释度/浓度应由使用者检定。

应 用	Ab评论	说明
IHC-P		1/100. Perform heat mediated antigen retrieval with citrate buffer pH 6 before commencing with IHC staining protocol.

靶标

功能

Probable oxidoreductase that has a dual role in controlling cellular life and death; during apoptosis, it is translocated from the mitochondria to the nucleus to function as a proapoptotic factor in a caspase-independent pathway, while in normal mitochondria, it functions as an antiapoptotic factor via its oxidoreductase activity. The soluble form (AlFsol) found in the nucleus induces 'parthanatos' i.e., caspase-independent fragmentation of chromosomal DNA. Interacts with ElF3G,and thereby inhibits the ElF3 machinery and protein synthesis, and activates casapse-7 to amplify apoptosis. Plays a critical role in caspase-independent, pyknotic cell death in hydrogen peroxide-exposed cells. Binds to DNA in a sequence-independent manner.

疾病相关

Defects in AIFM1 are the cause of combined oxidative phosphorylation deficiency type 6 (COXPD6) [MIM:300816]. It is a mitochondrial disease resulting in a neurodegenerative disorder characterized by psychomotor delay, hypotonia, areflexia, muscle weakness and wasting.

序列相似性

Belongs to the FAD-dependent oxidoreductase family.


翻译后修饰

Under normal conditions, a 54-residue N-terminal segment is first proteolytically removed during or just after translocation into the mitochondrial intermembrane space (IMS) by the mitochondrial processing peptidase (MPP) to form the inner-membrane-anchored mature form (AlFmit). During apoptosis, it is further proteolytically processed at amino-acid position 101 leading to the generation of the mature form, which is confined to the mitochondrial IMS in a soluble form (AlFsol). AlFsol is released to the cytoplasm in response to specific death signals, and translocated to the nucleus, where it induces nuclear apoptosis in a caspase-independent manner.

细胞定位

Mitochondrion intermembrane space. Mitochondrion inner membrane. Cytoplasm. Nucleus. Cytoplasm > perinuclear region. Proteolytic cleavage during or just after translocation into the mitochondrial intermembrane space (IMS) results in the formation of an inner-membrane-anchored mature form (AlFmit). During apoptosis, further proteolytic processing leads to a mature form, which is confined to the mitochondrial IMS in a soluble form (AlFsol). AlFsol is released to the cytoplasm in response to specific death signals, and translocated to the nucleus, where it induces nuclear apoptosis. Colocalizes with ElF3G in the nucleus and perinuclear region.

图片

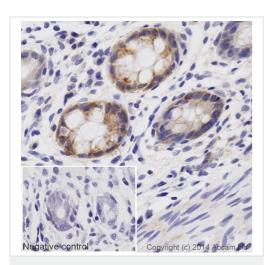


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Biotin Anti-AIF antibody [E20] - Mitochondrial Marker (ab197526)

IHC image of AIF staining in a section of formalin-fixed paraffinembedded Hu_Colon_Norm*, performed on a Leica Bond system using the standard protocol B. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab197526 at 1/100 dilution for 15 mins at room temperature and detected using an HRP conjugated ABC system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. The inset negative control image is taken from an identical assay without primary antibody.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

*Tissue obtained from the Human Research Tissue Bank, supported by the NIHR Cambridge Biomedical Research Centre

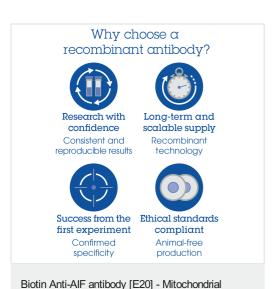

Marker (ab197526)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.cn/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors