# abcam

### Product datasheet

# Anti-Fibronectin antibody - BSA and Azide free ab23750

★★★★★ 19 Abreviews 195 References 2 图像

#### 概述

产**品名称** Anti-Fibronectin抗体- BSA and Azide free

描述 兔多克隆抗体to Fibronectin - BSA and Azide free

**宿主** Rabbit

 经测试应用
 适用于: IHC-P

 种属反应性
 与反应: Human

预测可用于: Mouse, Rat 4

免疫原 Full length native protein (purified) corresponding to Mouse Fibronectin.

常规说明 This antibody may be used to identify fibronectin during wound healing, during tumour

progression, and tumour invasion respectively.

The Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

#### 性能

形式 Liquic

**存放说明** Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.

存储溶液 Constituent: PBS

无载体 是

纯**度** Affinity purified

Primary antibody说明 This antibody may be used to identify fibronectin during wound healing, during tumour

progression, and tumour invasion respectively.

**克隆** 多克隆

**同种型** IgG

1

# The Abpromise guarantee

# Abpromise™承诺保证使用ab23750于以下的经测试应用

"应用说明"部分下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。


应 <b>用</b>	Ab评论	说明
IHC-P	★★★★☆ (6)	1/1000. Perform heat mediated antigen retrieval before commencing with IHC staining protocol.

靶标		
功能	Fibronectins bind cell surfaces and various compounds including collagen, fibrin, heparin, DNA, and actin. Fibronectins are involved in cell adhesion, cell motility, opsonization, wound healing, and maintenance of cell shape. Involved in osteoblast compaction through the fibronectin fibrillogenesis cell-mediated matrix assembly process, essential for osteoblast mineralization. Participates in the regulation of type I collagen deposition by osteoblasts.  Anastellin binds fibronectin and induces fibril formation. This fibronectin polymer, named superfibronectin, exhibits enhanced adhesive properties. Both anastellin and superfibronectin inhibit tumor growth, angiogenesis and metastasis. Anastellin activates p38 MAPK and inhibits lysophospholipid signaling.	
组 <b>织特异性</b>	Plasma FN (soluble dimeric form) is secreted by hepatocytes. Cellular FN (dimeric or cross-linked multimeric forms), made by fibroblasts, epithelial and other cell types, is deposited as fibrils in the extracellular matrix. Ugl-Y1, Ugl-Y2 and Ugl-Y3 are found in urine.	
疾病相关	Glomerulopathy with fibronectin deposits 2	
序列相似性	Contains 12 fibronectin type-I domains.  Contains 2 fibronectin type-II domains.  Contains 16 fibronectin type-III domains.	
发 <b>展</b> 阶段	Ugl-Y1, Ugl-Y2 and Ugl-Y3 are present in the urine from 0 to 17 years of age.	
翻译后修饰	Sulfated.  It is not known whether both or only one of Thr-2064 and Thr-2065 are/is glycosylated.  Forms covalent cross-links mediated by a transglutaminase, such as F13A or TGM2, between a glutamine and the epsilon-amino group of a lysine residue, forming homopolymers and heteropolymers (e.g. fibrinogen-fibronectin, collagen-fibronectin heteropolymers).  Phosphorylated by FAM20C in the extracellular medium.  Proteolytic processing produces the C-terminal NC1 peptide, anastellin.	

Secreted, extracellular space, extracellular matrix.

# 图片


细胞定位

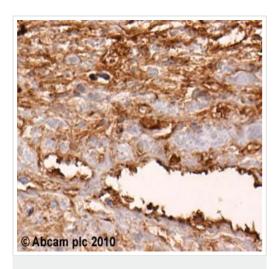

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Fibronectin antibody - BSA and Azide free (ab23750)

ab23750 (2ug/ml) staining Fibronectin in human aorta using an automated system (DAKO Autostainer Plus). Using this protocol there is strong staining Fibronectin in the cytoplasm and extracellularly in the connective tissue.

Sections were rehydrated and antigen retrieved with the Dako 3 in 1 AR buffer EDTA pH 9.0 in a DAKO PT link. Slides were peroxidase blocked in 3% H2O2 in methanol for 10 mins. They were then blocked with Dako Protein block for 10 minutes (containing casein 0.25% in PBS) then incubated with primary antibody for 20 min and detected with Dako envision flex amplification kit for 30 minutes. Colorimetric detection was completed with Diaminobenzidine for 5 minutes. Slides were counterstained with Haematoxylin and coverslipped under DePeX. Please note that, for manual staining, optimization of primary antibody concentration and incubation time is recommended. Signal amplification may be required.

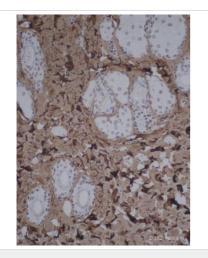

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Fibronectin antibody - BSA and Azide free (ab23750)

This image is courtesy of an anonymous abreview.

ab23750 staining Fibronectin in Rat skin tissue sections by Immunohistochemistry (IHC-P - paraformaldehyde-fixed, paraffinembedded sections). Tissue was fixed with formaldehyde and blocked with Immunoblock® (5 µg/ml) for 30 minutes at 25°C; antigen retrieval was by heat mediation in a citrate buffer. Samples were incubated with primary antibody (1/500 in diluent) for 1 hour 30 minutes at 25°C. An undiluted Goat anti-rabbit HRP/DAB polydetector polyclonal was used as the secondary antibody.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- · Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.cn/abpromise">https://www.abcam.cn/abpromise</a> or contact our technical team.

#### Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors