# abcam

## Product datasheet

# Anti-FOXP3 antibody [EPR15038-69] ab191416



重组 RabMAb

4 References 7 图像

概述

产品名称 Anti-FOXP3抗体[EPR15038-69]

描述 **兔**单**克隆抗体**[EPR15038-69] to FOXP3

宿主 Rabbit

适用于: IHC-Fr, IHC-P, WB 经测试应用

种属反应性 与反应: Human

免疫原 Recombinant fragment. This information is proprietary to Abcam and/or its suppliers.

阳性对照 IHC-P: FFPE Human normal tissue sections. Human Hodgkin lymphoma tissue. IHC-Fr: Human

tonsil tissue sections.

常规说明 This product is a recombinant monoclonal antibody, which offers several advantages including:

- High batch-to-batch consistency and reproducibility

- Improved sensitivity and specificity

- Long-term security of supply

- Animal-free production

For more information see here.

Our RabMAb® technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to **RabMAb**® **patents**.

性能

形式 Liquid

存放说明 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C long

term. Avoid freeze / thaw cycle.

存储溶液 pH: 7.2

Preservative: 0.01% Sodium azide

Constituents: 40% Glycerol, 59% PBS, 0.05% BSA

纯度 Protein A purified

克隆 单克隆

克隆编号 EPR15038-69

同种型 lgG

#### The Abpromise guarantee

#### Abpromise™承诺保证使用ab191416于以下的经测试应用

"应用说明"部分下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。

应用	Ab评论	说明
IHC-Fr		Use a concentration of 5 µg/ml.
IHC-P		1/100 - 1/250. Perform heat mediated antigen retrieval with Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.
WB		Use a concentration of 5 µg/ml. Detects a band of approximately 50 kDa (predicted molecular weight: 47 kDa).

#### 靶标

功能 Probable transcription factor. Plays a critical role in the control of immune response.

疾病相关 Defects in FOXP3 are the cause of immunodeficiency polyendocrinopathy, enteropathy, X-linked

syndrome (IPEX) [MIM:304790]; also known as X-linked autoimmunity-immunodeficiency syndrome. IPEX is characterized by neonatal onset insulin-dependent diabetes mellitus,


infections, secretory diarrhea, trombocytopenia, anemia and eczema. It is usually lethal in infancy.

序列相似性 Contains 1 C2H2-type zinc finger.

Contains 1 fork-head DNA-binding domain.

细胞定位 Nucleus.

### 图片

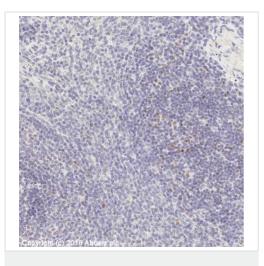


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-FOXP3 antibody
[EPR15038-69] (ab191416)

IHC image of FOXP3 staining in a section of formalin-fixed paraffinembedded normal human spleen\* performed on a Leica BOND<sup>TM</sup> system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20mins. The section was then incubated with ab191416, 5ug/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. The inset secondary-only control image is taken from an identical assay without primary antibody.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

<sup>\*</sup>Tissue obtained from the Human Research Tissue Bank,

supported by the NIHR Cambridge Biomedical Research Centre

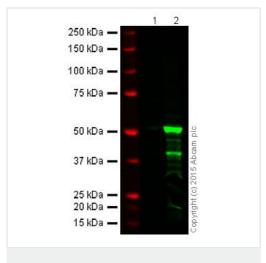

Western blot - Anti-FOXP3 antibody [EPR15038-69] (ab191416)

All lanes : Anti-FOXP3 antibody [EPR15038-69] (ab191416) at 5  $\,\mu g/ml$ 


Lane 1: HEK293T cell lysate

Lane 2: HEK293T cell lysate overexpressing Human FOXP3

Lysates/proteins at 20 µg per lane.

Performed under reducing conditions.

Predicted band size: 47 kDa

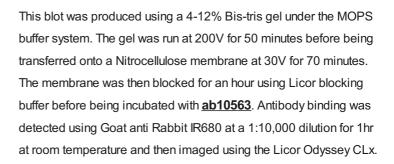




Immunohistochemistry (Frozen sections) - Anti-FOXP3 antibody [EPR15038-69] (ab191416)

IHC image of FOXP3 staining in a section of formalin-fixed frozen normal human tonsil\* performed on a Leica BOND<sup>TM</sup> system using the standard protocol F. The section was then incubated with ab191416, 5ug/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. The inset secondary-only control image is taken from an identical assay without primary antibody.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

\*Tissue obtained from the Human Research Tissue Bank, supported by the NIHR Cambridge Biomedical Research Centre

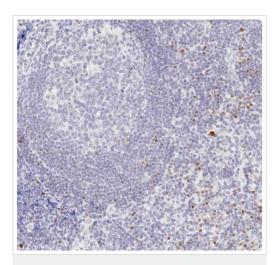


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-FOXP3 antibody
[EPR15038-69] (ab191416)

IHC image of FOXP3 staining in a section of formalin-fixed paraffinembedded normal human tonsil\* performed on a Leica BOND<sup>TM</sup> system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20mins. The section was then incubated with ab191416, 5ug/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. The inset secondary-only control image is taken from an identical assay without primary antibody.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.


\*Tissue obtained from the Human Research Tissue Bank, supported by the NIHR Cambridge Biomedical Research Centre

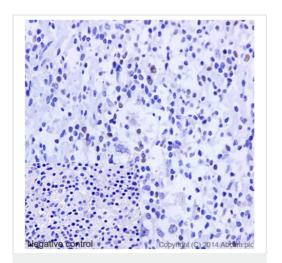

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-FOXP3 antibody
[EPR15038-69] (ab191416)


Immunohistochemical analysis of paraffin-embedded Human Hodgkin lymphoma tissue labeling FOXP3 with ab191416 at 1/250 dilution (0.9  $\mu$ g/ml) followed by pre-diluted HRP Polymer for Rabbit lgG secondary antibody and counter-stained with Hematoxylin. Inset: Negative control: using PBS instead of primary antibody.

Perform heat mediated antigen retrieval with Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

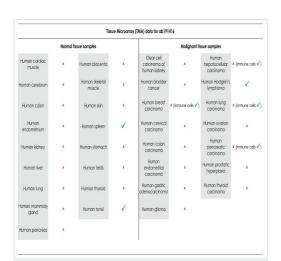

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-FOXP3 antibody
[EPR15038-69] (ab191416)

Tissue Microarrays stained for "Anti-FOXP3 antibody [EPR15038-69]" using "ab191416" in immunohistochemical analysis. This table provides a detailed overview of positive (tick mark) and negative (cross mark) staining per sample type tested. The sections were pre-treated using Heat mediated antigen retrieval using Bond™ Epitope Retrieval Solution 2 (pH 9.0) for 20 minutes. The sections were incubated with ab191416 for 30 mins at room temperature followed by a ready to use Rabbit specific IHC polymer detection kit HRP/DAB (ab209101). The immunostaining was performed on a Leica Biosystems BOND® RX instrument.

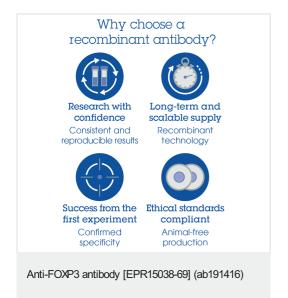

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

#### Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery

- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.cn/abpromise">https://www.abcam.cn/abpromise</a> or contact our technical team.

#### Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors