

Anti-NeuN antibody - Neuronal Marker ab104225

★★★★★ [11 Abreviews](#) [224 References](#) [9 图像](#)

概述	
产品名称	Anti-NeuN抗体- Neuronal Marker
描述	兔多克隆抗体to NeuN - Neuronal Marker
宿主	Rabbit
经测试应用	适用于: WB, IHC-FrFI, IHC-P, IHC-FoFr
种属反应性	与反应: Mouse, Rat, Human
免疫原	Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.
常规说明	<p>The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.</p> <p>If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&As</p>
性能	
形式	Liquid
存放说明	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C. Avoid freeze / thaw cycle.
存储溶液	Preservative: 0.03% Sodium azide Constituents: 49.98% PBS, 49.98% Glycerol (glycerin, glycerine)
纯度	Affinity purified
克隆	多克隆
同种型	IgG
应用	

The Abpromise guarantee [Abpromise™](#)承诺保证使用ab104225于以下的经测试应用


“应用说明”部分 下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。

应用	Ab评论	说明
WB		1/5000 - 1/10000. Predicted molecular weight: 48 kDa.
IHC-FrFI	★★★★★ (1)	Use at an assay dependent concentration.
IHC-P	★★★★★★ (6)	1/500 - 1/1000.
IHC-FoFr	★★★★★ (1)	Use at an assay dependent concentration.

靶标


功能	RNA-binding protein that regulates alternative splicing events.
序列相似性	Contains 1 RRM (RNA recognition motif) domain.
细胞定位	Nucleus. Cytoplasm.

图片

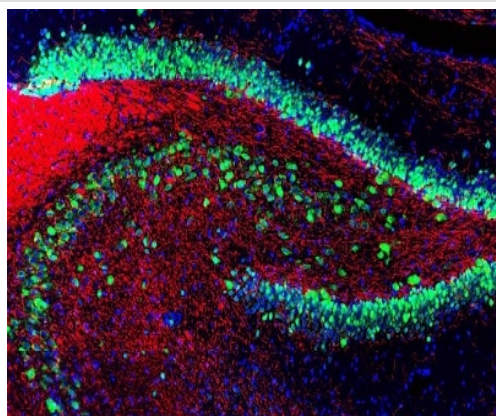

Immunohistochemical analysis of paraformaldehyde fixed frozen section Rat hippocampus tissue labeling NeuN with ab104225 at 1/1000 dilution.


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-NeuN antibody - Neuronal Marker (ab104225)

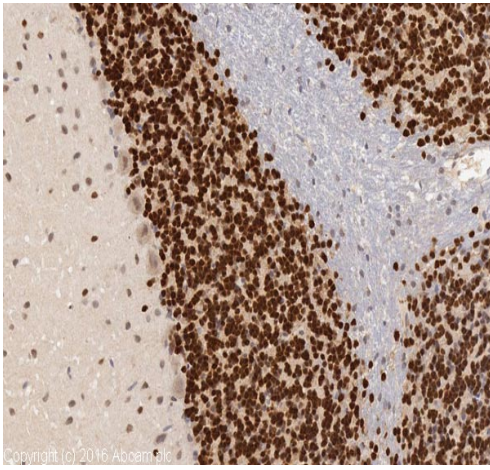

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-NeuN antibody - Neuronal Marker (ab104225)

IHC image of NeuN staining in rat cerebellum formalin fixed paraffin embedded tissue section, performed on a Leica Bond™ system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab104225, 1in500 dilution, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.


For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

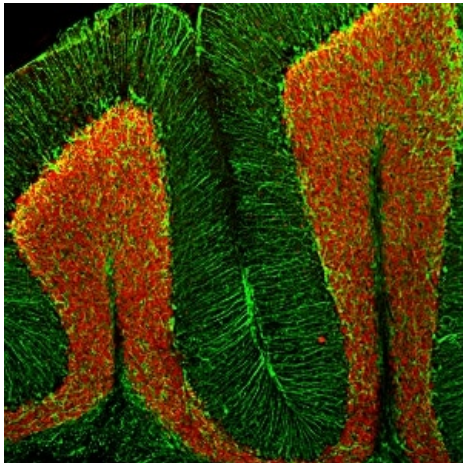

Immunohistochemistry - Free Floating - Anti-NeuN antibody - Neuronal Marker (ab104225)

Immunofluorescent analysis of a section of adult mouse cerebellum stained with ab104225 at a dilution 1:5,000 in red, co-stained with chicken pAb to GFAP, dilution 1:5,000, in green.

Following transcardial perfusion of mouse with 4% paraformaldehyde, brain was post fixed for 24 hours, cut to 45 µM, and free-floating sections. ab104225 stains the nuclei of neurons in the cerebellar granule layer. The GFAP antibody stains the processes of Bergmann glia in the molecular layer and astroglia in the granule and white matter layers.

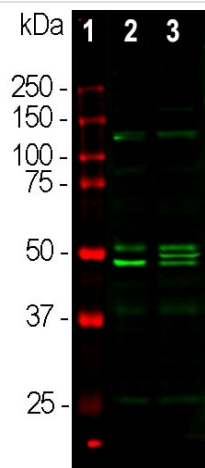

Western blot - Anti-NeuN antibody - Neuronal Marker (ab104225)

All lanes : Anti-NeuN antibody - Neuronal Marker (ab104225) at 1/1000 dilution

Lane 2 : Mouse whole brain lysate

Lane 3 : Rat whole brain lysate


Predicted band size: 48 kDa



Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-NeuN antibody - Neuronal Marker (ab104225)

IHC image of NeuN staining in human cerebellum formalin fixed paraffin embedded tissue section, performed on a Leica Bond™ system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab104225, 1in500 dilution, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.


For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

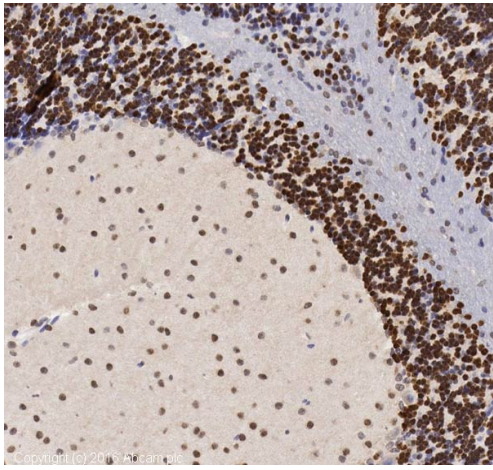

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-NeuN antibody - Neuronal Marker (ab104225)

IHC image of NeuN staining in mouse cerebellum formalin fixed paraffin embedded tissue section, performed on a Leica Bond™ system using the standard protocol B. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab104225, 1in500 dilution, for 15 mins at room temperature. A goat anti-rabbit biotinylated secondary antibody was used to detect the primary, and visualized using an HRP conjugated ABC system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.


For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

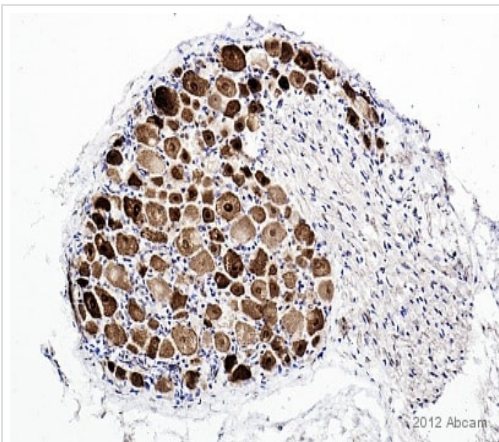

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-NeuN antibody - Neuronal Marker (ab104225)


This image is courtesy of an abreview submitted by Carl Hobbs, King's College London, United Kingdom

IHC-P image of FOX3/NeuN staining using rat DRG using **ab104225** (1:500). The sections were subjected to heat mediated antigen retrieval using citric acid pH 6. The sections were blocked using 1% BSA for 10 mins at 21°C. **ab104225** was incubated for 2 hours at 21°C. The secondary antibody used was Goat polyclonal to Rabbit IgG conjugated to biotin (1:200).

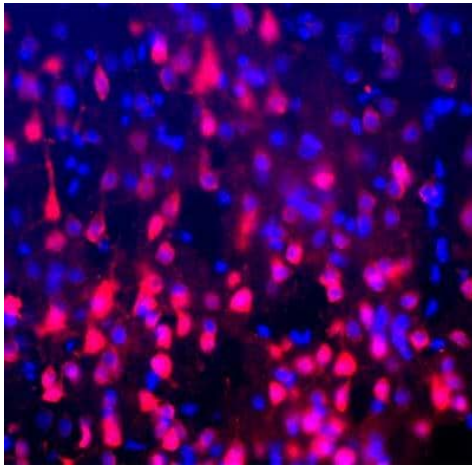

Immunohistochemistry (PFA perfusion fixed frozen sections) - Anti-NeuN antibody - Neuronal Marker (ab104225)

ab104225 at 1/1000 dilution, staining FOX3 (red) in rat brain cortex and counterstained for DNA in blue.

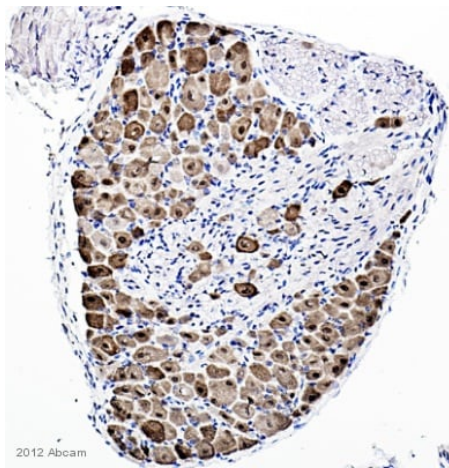

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-NeuN antibody - Neuronal Marker (ab104225)

This image is courtesy of an abreview submitted by Carl Hobbs, King's College London, United Kingdom

IHC-P image of FOX3/NeuN staining using mouse DRG using **ab104225** (1:500). The sections were subjected to heat mediated antigen retrieval using citric acid pH 6. The sections were blocked using 1% BSA for 10 mins at 21°C. **ab104225** was incubated for 2 hours at 21°C. The secondary antibody used was Goat polyclonal to Rabbit IgG conjugated to biotin (1:200).

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.cn/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors