


Anti-Periostin antibody ab79946

★★★★★ [2 Abreviews](#) [10 References](#) [3 图像](#)

概述

产品名称	Anti-Periostin抗体
描述	兔多克隆抗体to Periostin
宿主	Rabbit
经测试应用	适用于: IHC-P, WB, ICC/IF
种属反应性	与反应: Human 预测可用于: Mouse, Chicken, Cow 
免疫原	Synthetic peptide corresponding to Human Periostin aa 250-350 conjugated to keyhole limpet haemocyanin. (Peptide available as ab90847)
阳性对照	This antibody gave a positive signal in WB in Human breast tumour tissue lysate. This antibody gave a positive result in IHC in the following FFPE tissue: Human breast fibroadenocarcinoma.
常规说明	<p>The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets your needs before purchasing.</p> <p>If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be found below, along with publications, customer reviews and Q&As</p>

性能

形式	Liquid
存放说明	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C or -80°C. Avoid freeze / thaw cycle.
存储溶液	pH: 7.40 Preservative: 0.02% Sodium azide Constituent: PBS
	Batches of this product that have a concentration < 1mg/ml may have BSA added as a stabilising agent. If you would like information about the formulation of a specific lot, please contact our scientific support team who will be happy to help.

纯度	Immunogen affinity purified
克隆	多克隆
同种型	IgG

应用

The Abpromise guarantee **Abpromise™**承诺保证使用ab79946于以下的经测试应用


“应用说明”部分 下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。

应用	Ab评论	说明
IHC-P	★★★★★ (2)	Use a concentration of 0.1 µg/ml.
WB		Use a concentration of 1 µg/ml. Detects a band of approximately 80, 86, 86, 93 kDa (predicted molecular weight: 93 kDa). This antibody was raised against an immunogen that is predicted to recognize isoforms 1,2,3 & 4 of Human Periostin.
ICC/IF		Use a concentration of 1 µg/ml.

靶标

功能	Binds to heparin. Induces cell attachment and spreading and plays a role in cell adhesion. May play a role in extracellular matrix mineralization.
组织特异性	Widely expressed with highest levels in aorta, stomach, lower gastrointestinal tract, placenta, uterus and breast. Up-regulated in epithelial ovarian tumors. Not expressed in normal ovaries. Also highly expressed at the tumor periphery of lung carcinoma tissue but not within the tumor. Overexpressed in breast cancers.
序列相似性	Contains 1 EMI domain. Contains 4 FAS1 domains.
翻译后修饰	Gamma-carboxyglutamate residues are formed by vitamin K dependent carboxylation. These residues are essential for the binding of calcium.
细胞定位	Secreted > extracellular space > extracellular matrix.


图片

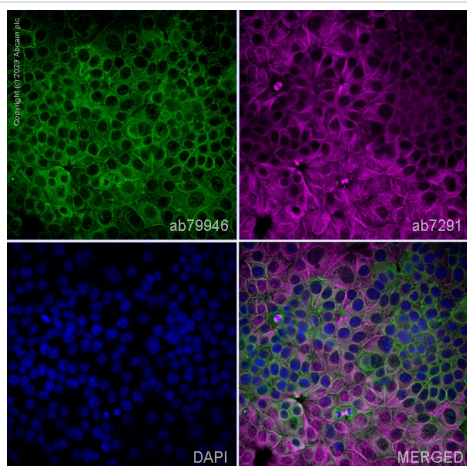

Immunocytochemistry/ Immunofluorescence - Anti-Periostin antibody (ab79946)

ab79946 staining Periostin in MCF7 cells. The cells were fixed with 100% methanol (5 min), permeabilized with 0.1% PBS-Tween for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1%PBS-Tween for 1h. The cells were then incubated overnight at 4°C with ab79946 at 5µg/ml and **ab7291**, Mouse monoclonal [DM1A] to alpha Tubulin - Loading Control. Cells were then incubated with **ab150081**, Goat polyclonal Secondary Antibody to Rabbit IgG - H&L (Alexa Fluor® 488), pre-adsorbed at 1/1000 dilution (shown in green) and **ab150120**, Goat polyclonal Secondary Antibody to Mouse IgG - H&L (Alexa Fluor® 594), pre-adsorbed at 1/1000 dilution (shown in pseudocolour magenta). Nuclear DNA was labelled with DAPI (shown in blue).


Image was acquired with a high-content analyser (Operetta CLS, Perkin Elmer) and a maximum intensity projection of confocal sections is shown.

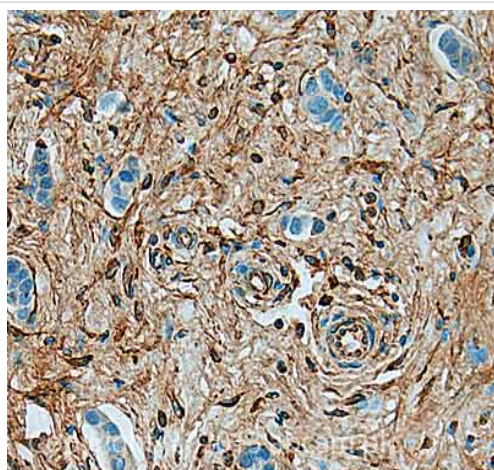

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Periostin antibody (ab79946)

IHC image of Periostin staining in Human breast fibroadenocarcinoma formalin fixed paraffin embedded tissue section, performed on a Leica Bond™ system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab79946, 0.1µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

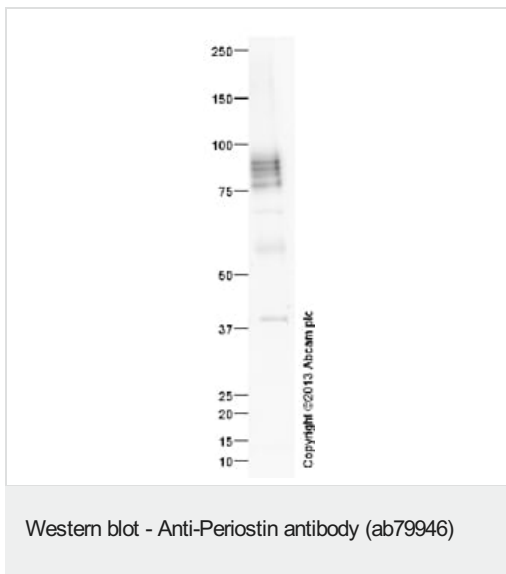

Anti-Periostin antibody (ab79946) at 1 µg/ml (blocked with 3% milk)
+ Breast (Human) Cytoplasmic Lysate - tumor tissue (**ab30085**) at
10 µg

Secondary

Goat polyclonal to Rabbit IgG - H&L - Pre-Adsorbed (HRP) at
1/10000 dilution

Developed using the ECL technique.

Performed under reducing conditions.

Predicted band size: 93 kDa

Observed band size: 80,80,86,87,93 kDa

Additional bands at: 38 kDa. We are unsure as to the identity of
these extra bands.

Exposure time: 20 minutes

This antibody was raised against an immunogen that is predicted to
recognize isoforms 1,2,3 & 4 of Human Periostin. The predicted
molecular weights of isoforms 1,2,3 & 4 are 93, kDa, 87 kDa, 87
kDa and 83 kDa respectively.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.cn/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors