

Human Angiogenesis Antibody Array B (30 Targets)- Quantitative ab197421

3 图像

概述

产品名称	人Angiogenesis抗体阵列B (30 Targets)- Quantitative
样品类型	Serum, Plasma, Other biological fluids, Cell Lysate, Cell culture media, Tissue Lysate
检测类型	Quantitative
检测时间	6h 00m
种属反应性	与反应: Human
产品概述	ab197421 can be used for quantitative measurement of 30 Human Angiogenesis proteins. Suitable for serum, plasma, cell culture media, other body fluids, cell and tissue lysates. Targets: ANG-1, Angiostatin, CXCL16, EGF, FGF-4, Follistatin, G-CSF, GM-CSF, I-309, IL-1 β , IL-4, IL-10, IL-12p40, IL-12p70, ITAC, MCP-2, MCP-3, MCP-4, MMP-1, MMP-9, PECAM-1, TGF α , TGF β 3, Tie-1, Tie-2, uPAR, VEGF, VEGF R2, VEGF R3, VEGF-D.

[Learn more about Quantitative antibody arrays](#)

说明

Quantitative antibody arrays can be used to quantitate up to 40 cytokines with as little as 50 μ L of sample. Arrays are available for 400 human and 200 mouse proteins.

Each glass slide is spotted with 16 identical antibody arrays and is provided with a 16 well gasket to allow separate samples to be applied to each array. Each antibody is spotted in quadruplicate in each array. For reproducible quantitation, eight of the arrays are used with a cocktail of protein standards to produce a standard curve. The same 8 standard curve arrays can be used across multiple slides, allowing measurement of 22 experimental samples with 2 slides, 36 with 3 slides etc.

For high through-put, 4 slides can be nested into a tray matching a standard microplate, allowing automated processing with a liquid handling workstation.

Array processing can be completed within one working day (see workflow diagram below) and arrays can be analyzed with a wide number of laser glass slide array / gene microarray scanners. If you don't have a suitable scanner then we recommend our membrane antibody arrays, which can be analyzed with any WB chemiluminescent reader.

The standard is affinity purified.

[Learn more about quantitative antibody arrays](#)

[Learn more about membrane antibody arrays](#)


性能

存放说明

Store at -20°C. Please refer to protocols.

组件	50 tests	22 tests	8 tests
20X Wash Buffer I	6 x 30ml	3 x 30ml	2 x 30ml
20X Wash Buffer II	2 x 30ml	1 x 30ml	1 x 30ml
Cy3 equivalent dye-conjugated Streptavidin	4 x 5µl	2 x 5µl	1 x 5µl
Human Angiogenesis Array B (30 T) Biotinylated Antibody Cocktail	4 units	2 units	1 unit
Human Angiogenesis Array B (30 T) Glass Slide	4 slides	2 slides	1 slide
Human Angiogenesis Array B (30 T) Lyophilized Cytokine Standard Mix	2 vials	2 vials	1 vial
Sample Diluent	2 x 15ml	1 x 15ml	1 x 15ml
Slide Washer/Dryer (30 mL Centrifuge Tube)	2 units	1 unit	1 unit

图片

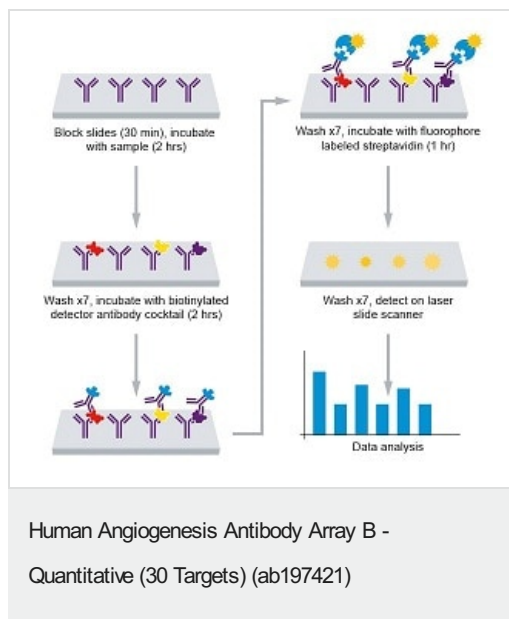


Assay summary for Abcam Human Angiogenesis Antibody Array B (30 Targets) - Quantitative (ab197421)

POS1	POS2	NEG
ANG-1	Angiostatin	CXCL16
EGF	FGF-4	Follistatin
G-CSF	GM-CSF	I-309
IL-1 β	IL-4	IL-10
IL-12p40	IL-12p70	I-TAC
MCP-2	MCP-3	MCP-4
MMP-1	MMP-9	PECAM-1
TGF α	TGF β 3	Tie-1
Tie-2	uPAR	VEGF
VEGF R2	VEGF R3	VEGF-D

Array Map for Abcam Human Angiogenesis Antibody Array B (30 Targets) - Quantitative (ab197421)

Human Angiogenesis Antibody Array B -
Quantitative (30 Targets) (ab197421)

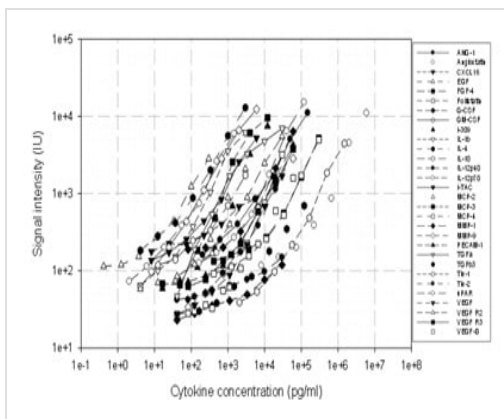

Standard Curves obtained from Abcam Human Angiogenesis Antibody Array B (30 Targets) - Quantitative (ab197421)

Human Angiogenesis Antibody Array B -
Quantitative (30 Targets) (ab197421)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.cn/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors