

# Mounting Medium With DAPI - Aqueous, Fluoroshield ab104139

★★★★★ [8 Abreviews](#) [327 References](#) [2 图像](#)

### 概述

#### 产品名称

封固剂With DAPI - Aqueous, Fluoroshield

#### 经测试应用

**适用于:** IHC-Fr, ICC/IF, IHC-P, In situ hybridization

#### 常规说明

Mounting Medium with DAPI ab104139 (previously called Fluoroshield Mounting Medium With DAPI) is an aqueous anti-fade mounting medium designed to preserve fluorescence when imaging tissues and cell samples. It is fortified with DAPI as a counter-stain for DNA (DAPI also stains RNA).

The formulation prevents rapid photobleaching of FITC, Texas Red, AMCA, Cy2, Cy3, Cy5, Alexa fluoro 488, Alexa fluoro 594, Green fluorescent protein (GFP), tetramethyl rhodamine, R-Phycoerythrin (R-PE), Phyocyanin (PC), and Allophycocyanin (APC).

Fluorescence is retained during prolonged storage at 4°C in the dark. This medium does not contain phenylenediamine, which destroys immunofluorescence of Cy dyes, R-PE, PC and APC.

#### Mounting media products

For fluorescent cell and tissue staining, Abcam recommends aqueous, anti-fade **Fluorescence Mounting Medium ab104135**, **Mounting Medium with PI ab104129**, and this product, Mounting Medium with DAPI ab104139. For thick sections or tissues containing lots of fat, we recommend **Glycerol Mounting Medium with DAPI ab188804**.

For chromogenic immunohistochemistry, such as with DAB, AEC, or Fast Red, we recommend aqueous **Mounting Medium ab64230**, or organic **Limonene Mounting Medium ab104141**.

#### Procedure

- Bring the vial to room temperature.
- Rinse slide to be mounted with distilled or deionized water, touch the edges of slide on a paper towel to remove excess water. Place slides on a flat surface.
- Turn the vial upside down and open the dropper to remove any air bubbles.
- Apply 3-4 drops of mounting medium directly on top of the specimen and spread out evenly by tilting slide back and forth or spread evenly with a 0.2 ml plastic pipette tip making sure the tissue is not touched. Excess medium can be removed by touching the edges of slide against paper towel.
- Let stand at room temperature for about 5 minutes.
- Apply cover slip carefully avoiding air bubbles.

- The specimen is ready for visualization under a microscope.

Seal the edges of cover slip with nail polish, any organic medium or **Limonene mounting medium ab104141**. If a coverslip is not sealed, air bubbles will appear in a few days.

Store the slide in the dark at 2-8°C.

#### Removal of Coverslip

Coverslip can be removed before sealing the edges. Soak slide in warm (37°C) distilled or deionized water for several minutes. Carefully and slowly move the coverslip. Soak in water for an additional few minutes to remove coverslip. Rinse slide several times with warm water to remove all mounting medium. The slide can be remounted again.

#### Refractive Index

1.364 ± 0.002. This number applies to this mounting medium in solution. Refractive indexes change when the water solvent evaporates and mounting media dries on slides. We do not have the means to measure the refractive indexes of dry mounting mediums; however we expect the numbers to go higher when dried. The refractive index of water is 1.3330.

### 性能

形式	Liquid
存放说明	Store at +4°C. Do Not Freeze. Store In the Dark.
存储溶液	Preservative: 0.09% Sodium azide Constituents: 0.0002% 4',6-Diamidino-2-phenylindole dihydrochloride, 84% Water


### 应用

The Abpromise guarantee **Abpromise™**承诺保证使用ab104139于以下的经测试应用

“应用说明”部分 下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。

应用	Ab评论	说明
IHC-Fr		Use at an assay dependent concentration.
ICC/IF		Use at an assay dependent concentration.
IHC-P		Use at an assay dependent concentration.
In situ hybridization		Use at an assay dependent concentration.


### 图片

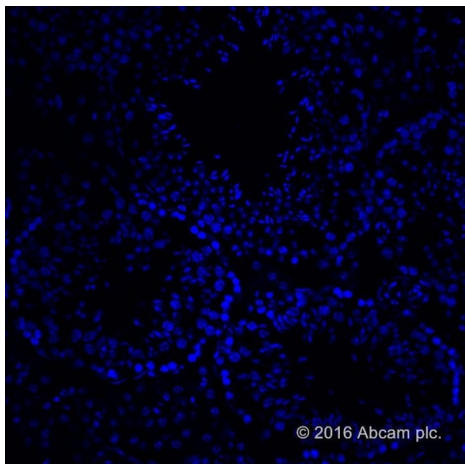

Immunohistochemistry (Frozen sections) -  
Fluoroshield Mounting Medium With DAPI  
(ab104139)

This image is courtesy of an abreview submitted by  
Bryan Niedenberger.

Adult mouse testes were fixed with 4% paraformaldehyde. Tissue was embedded in O.C.T. and frozen. 5 micron sections were cut and transferred to slides. Sections were permeabilized with 0.1% Triton X-100 in PBS and stained with undiluted Fluoroshield Mounting Medium with DAPI (ab104139).

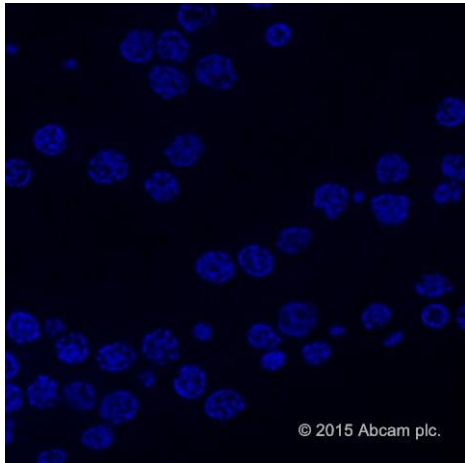

Immunocytochemistry/ Immunofluorescence -  
Fluoroshield Mounting Medium With DAPI  
(ab104139)

This image is courtesy of an Abreview submitted by  
Sebastian Alvarado.

ab104139 worked exactly as it should in resin sections and stained heterochromatic areas of the nucleus.

**Please note:** All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

**Our Abpromise to you: Quality guaranteed and expert technical support**

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.cn/abpromise> or contact our technical team.

#### **Terms and conditions**

---

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors