abcam

Product datasheet

FITC Anti-EpCAM antibody [B29.1 (VU-1D9)] ab8666

★★★★★ 2 Abreviews 19 References 1 图像

概述

产品名称 FITC荧光Anti-EpCAM抗体[B29.1 (VU-1D9)]

描述 FITC荧光小鼠单克隆抗体[B29.1 (VU-1D9)] to EpCAM

宿主 Mouse

偶联物 FITC. Ex: 493nm, Em: 528nm

特异性 This antibody reacts with a 40 kD cell surface glycoprotein called ESA.

经测试应用 适用于: ICC, IHC-P, WB, IP, ICC/IF

种属反应性 与反应: Human

免疫原 Tissue, cells or virus corresponding to Human EpCAM. Raised against carcinoma cell line of

human origin.

Database link: P16422

阳性对照 Breast, colon carcinoma and tonsil

常规说明 Dilute according to the particular application being used. In general, 0.05M borate pH 8.0

containing 0.15M sodium chloride, 0.02% sodium azide, is a good dilutent to use with most antibodies. Avoid diluting the entire contents of the vial at once since the diluted solution may have

reduced stability.

Fluorescein (FITC) conjugated at a ratio of 6 moles of Fluorescein to 1 mole of antibody.

Maximum excitation of FITC occurs at wavelength 492 nm and maximum emission occurs at

wavelength 520 nm.

The Life Science industry has been in the grips of a reproducibility crisis for a number of years. Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

性能

形式 Liquid

存放说明 Shipped at 4°C. Store at +4°C. Please see notes section.

存储溶液 pH: 7.40

1

Preservative: 0.1% Sodium azide

Constituent: 0.01% PBS

纯**度** Protein G purified

克隆 单克隆

克隆编号 B29.1 (VU-1D9)

同种型 lgG1 轻链类型 kappa

应用

The Abpromise guarantee Abpromise™承诺保证使用ab8666于以下的经测试应用

"应用说明"部分下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。

应用	Ab评论	说明
ICC		Use at an assay dependent concentration. PubMed: 17417652
IHC-P		Use at an assay dependent concentration. PubMed: 20421921
WB		Use at an assay dependent concentration. Predicted molecular weight: 35 kDa. This product was not quality controlled in Flow Cytometry applications. However, it may be applicable for Flow Cytometry at a concentration of 1-2 µg in 0.1 ml containing 10^6 cells.
IP		Use at an assay dependent concentration.
ICC/IF	★★★★☆ (1)	Use a concentration of 5 - 10 μg/ml. See Abreview.

靶标

功能 May act as a physical homophilic interaction molecule between intestinal epithelial cells (IECs)

and intraepithelial lymphocytes (IELs) at the mucosal epithelium for providing immunological barrier as a first line of defense against mucosal infection. Plays a role in embryonic stem cells proliferation and differentiation. Up-regulates the expression of FABP5, MYC and cyclins A and E.

组织特异性 Highly and selectively expressed by undifferentiated rather than differentiated embryonic stem

cells (ESC). Levels rapidly diminish as soon as ESC's differentiate (at protein levels). Expressed in almost all epithelial cell membranes but not on mesodermal or neural cell membranes. Found

on the surface of adenocarcinoma.

疾病相关 Defects in EPCAM are the cause of diarrhea type 5 (DIAR5) [MIM:613217]. It is an intractable diarrhea of infancy characterized by villous atrophy and absence of inflammation, with intestinal

epithelial cell dysplasia manifesting as focal epithelial tufts in the duodenum and jejunum. Defects in EPCAM are a cause of hereditary non-polyposis colorectal cancer type 8 (HNPCC8) [MIM:613244]. HNPCC is a disease associated with marked increase in cancer susceptibility. It is characterized by a familial predisposition to early-onset colorectal carcinoma (CRC) and extracolonic tumors of the gastrointestinal, urological and female reproductive tracts. HNPCC is reported to be the most common form of inherited colorectal cancer in the Western world.

2

沙开

Clinically, HNPCC is often divided into two subgroups. Type I is characterized by hereditary predisposition to colorectal cancer, a young age of onset, and carcinoma observed in the proximal colon. Type II is characterized by increased risk for cancers in certain tissues such as the uterus, ovary, breast, stomach, small intestine, skin, and larynx in addition to the colon. Diagnosis of classical HNPCC is based on the Amsterdam criteria: 3 or more relatives affected by colorectal cancer, one a first degree relative of the other two; 2 or more generation affected; 1 or more colorectal cancers presenting before 50 years of age; exclusion of hereditary polyposis syndromes. The term 'suspected HNPCC' or 'incomplete HNPCC' can be used to describe families who do not or only partially fulfill the Amsterdam criteria, but in whom a genetic basis for colon cancer is strongly suspected. Note=HNPCC8 results from heterozygous deletion of 3-prime exons of EPCAM and intergenic regions directly upstream of MSH2, resulting in transcriptional read-through and epigenetic silencing of MSH2 in tissues expressing EPCAM.

序列相似性

Belongs to the EPCAM family.

Contains 1 thyroglobulin type-1 domain.

翻译后修饰

Hyperglycosylated in carcinoma tissue as compared with autologous normal epithelia.


Glycosylation at Asn-198 is crucial for protein stability.

细胞定位

 $\label{lem:lembrane} Lateral\ cell\ membrane.\ Cell\ junction > tight\ junction.\ Co-localizes\ with\ CLDN7\ at\ the\ lateral\ cell$

membrane and tight junction.

图片

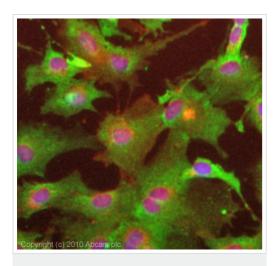

Immunocytochemistry/ Immunofluorescence - FITC Anti-EpCAM antibody [B29.1 (VU-1D9)] (ab8666)

ICC/IF image of ab8666 stained HepG2 cells. The cells were 4% formaldehyde fixed (10 min) and then incubated in 1%BSA / 10% normal goat serum / 0.3M glycine in 0.1% PBS-Tween for 1h to permeabilise the cells and block non-specific protein-protein interactions. The cells were then incubated with the antibody (ab8666, 5 μ g/ml) overnight at +4°C. Alexa Fluor® 594 WGA was used to label plasma membranes (red) at a 1/200 dilution for 1h. DAPI was used to stain the cell nuclei (blue) at a concentration of 1.43 μ M.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish

- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.cn/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors