

Anti-Granzyme B antibody [EPR20129-217] ab208586

重组 RabMAb

3 References **10 图像**

概述

产品名称	Anti-Granzyme B抗体[EPR20129-217]
描述	兔单克隆抗体[EPR20129-217] to Granzyme B
宿主	Rabbit
经测试应用	适用于: WB, IHC-P, mIHC 不适用于: Flow Cyt
种属反应性	与反应: Human
免疫原	Recombinant full length protein. This information is proprietary to Abcam and/or its suppliers.
阳性对照	WB: Human Granzyme B recombinant protein; IHC-P: Human colon and cervix cancer tissues. mIHC: Human breast cancer tissue.
常规说明	This product is a recombinant monoclonal antibody, which offers several advantages including: <ul style="list-style-type: none"> - High batch-to-batch consistency and reproducibility - Improved sensitivity and specificity - Long-term security of supply - Animal-free production For more information see here . Our RabMAb [®] technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to RabMAb[®] patents .

性能

形式	Liquid
存放说明	Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C long term. Avoid freeze / thaw cycle.
存储溶液	pH: 7.2 Preservative: 0.01% Sodium azide Constituents: 59% PBS, 40% Glycerol (glycerin, glycerine), 0.05% BSA
纯度	Protein A purified
克隆	单克隆
克隆编号	EPR20129-217

同种型IgG

应用

The Abpromise guarantee Abpromise™承诺保证使用ab208586于以下的经测试应用
“应用说明”部分 下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。


应用	Ab评论	说明
WB		1/10000. Detects a band of approximately 28 kDa (predicted molecular weight: 28 kDa).
IHC-P		1/250. Perform heat mediated antigen retrieval with Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.
mlHC		Use at an assay dependent concentration.

应用说明Is unsuitable for Flow Cyt.

靶标

功能	This enzyme is necessary for target cell lysis in cell-mediated immune responses. It cleaves after Asp. Seems to be linked to an activation cascade of caspases (aspartate-specific cysteine proteases) responsible for apoptosis execution. Cleaves caspase-3, -7, -9 and 10 to give rise to active enzymes mediating apoptosis.
序列相似性	Belongs to the peptidase S1 family. Granzyme subfamily. Contains 1 peptidase S1 domain.
细胞定位	Cytoplasmic granule. Cytoplasmic granules of cytolytic T-lymphocytes and natural killer cells.

图片

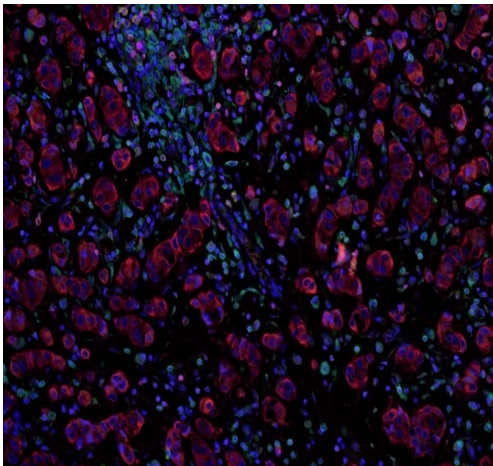

Multiplex immunohistochemistry - Anti-Granzyme B antibody [EPR20129-217] (ab208586)

This image is courtesy of ImmunoAtlas.

Fluorescence multiplex immunohistochemical analysis of Human breast cancer tissue (formalin-fixed paraffin-embedded section).

Merged staining of Anti-PD-L1 ([ab251611](#); cyan; Opal™ 520), Anti-Granzyme B ([ab219803](#); yellow; Opal™ 540), Anti-PD1 ([ab251613](#); magenta; Opal™ 570), Anti-pan Cytokeratin ([ab264485](#); red; Opal™ 620), Anti-EpCAM ([ab225894](#); red; Opal™ 620), Anti-CD8 alpha ([ab251596](#); green; Opal™ 650) and Anti-FOXP3 ([ab96048](#); orange; Opal™ 690). EpCAM and pan-cytokeratin share the same dye and color. Dyes are pseudo-colored for better contrast of the markers.


The immunostaining was performed on a Leica Biosystems BOND® MAX instrument with an Opal™ 6-Plex Detection Kit (NEL821001KT, Akoya Biosciences®).

The section was incubated in six rounds of staining; sequentially for

ab251611 (1/750 dilution), **ab219803** (1/250 dilution), **ab251613** (1/750 dilution), **ab264485** (0.5 µg/ml), **ab225894** (1/1250 dilution), **ab251596** (1/1500 dilution) and **ab96048** (10 µg/ml); each using a separate fluorescent tyramide signal amplification system. EDTA based antigen retrieval (Leica Biosystems BOND® Epitope Retrieval Solution 2, pH 9.0, 20 minutes) was used in between rounds of tyramide signal amplification to remove the antibody from the previous round, to avoid any cross-reactivity. DAPI (dark blue) was used as a nuclear counter stain.

Microscopy and pseudocoloring of individual Opal™ dyes was performed using a Vectra 3 Imaging System (Akoya Biosciences®).

This data is courtesy of ImmunoAtlas and it can be found [here](#).

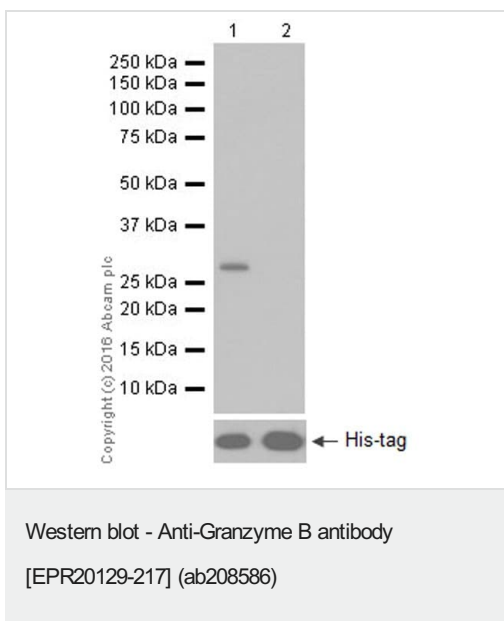

All lanes : Anti-Granzyme B antibody [EPR20129-217] (ab208586) at 1/10000 dilution

Lane 1 : Human Granzyme B recombinant protein

Lane 2 : Human Granzyme H recombinant protein

Lysates/proteins at 0.01 µg per lane.

Secondary

All lanes : Goat Anti-Rabbit IgG H&L (HRP) (**ab97051**) at 1/100000 dilution


Predicted band size: 28 kDa

Observed band size: 26 kDa

Exposure time: 1 second

Blocking/Dilution buffer: 5% NFDM/TBST.

Human Granzyme B and human Granzyme H recombinant protein contain aa21-247 and aa21-246 with a His-tag. These two recombinant proteins were made in house.

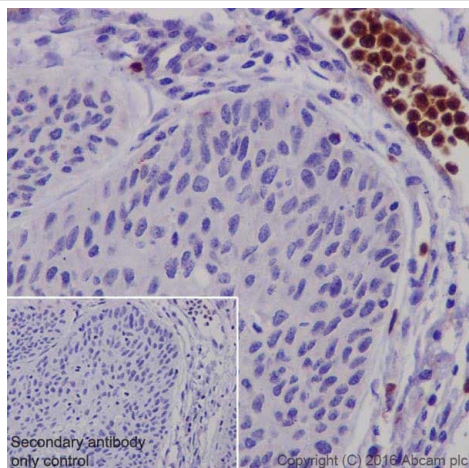


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Granzyme B antibody [EPR20129-217] (ab208586)

Immunohistochemical analysis of paraffin-embedded human cervix cancer tissue labeling Granzyme B with ab208586 at 1/250 dilution, followed by Goat Anti-Rabbit IgG H&L (HRP) ([ab97051](#)) at 1/500 dilution. Cytoplasmic staining on neutrophils and stroma cells of human cervix cancer is observed [PMID: 14512315]. Counter stained with Hematoxylin.

Secondary antibody only control: Used PBS instead of primary antibody, secondary antibody is Goat Anti-Rabbit IgG H&L (HRP) ([ab97051](#)) at 1/500 dilution.

Perform heat mediated antigen retrieval with Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

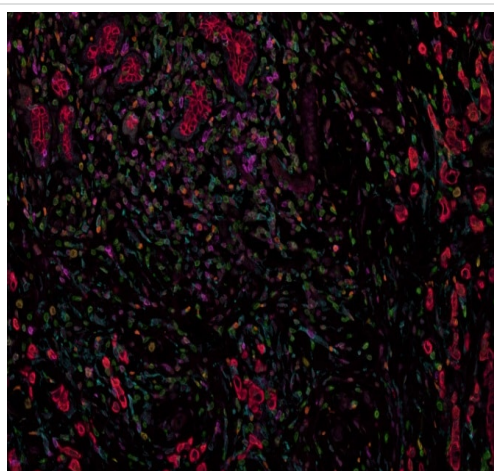

Multiplex immunohistochemistry - Anti-Granzyme B antibody [EPR20129-217] (ab208586)
This image is courtesy of ImmunoAtlas.

Fluorescence multiplex immunohistochemical analysis of Human breast cancer tissue (formalin-fixed paraffin-embedded section).

Merged staining of Anti-PD-L1 ([ab251611](#); cyan; Opal™ 520), Anti-Granzyme B ([ab219803](#); yellow; Opal™ 540), Anti-PD1 ([ab251613](#); magenta; Opal™ 570), Anti-pan Cytokeratin ([ab264485](#); red; Opal™ 620), Anti-EpCAM ([ab225894](#); red; Opal™ 620), Anti-CD8 alpha ([ab251596](#); green; Opal™ 650) and Anti-FOXP3 ([ab96048](#); orange; Opal™ 690). EpCAM and pan-cytokeratin share the same dye and color. Dyes are pseudo-colored for better contrast of the markers.


The immunostaining was performed on a Leica Biosystems BOND® MAX instrument with an Opal™ 6-Plex Detection Kit (NEL821001KT, Akoya Biosciences®).

The section was incubated in six rounds of staining; sequentially for [ab251611](#) (1/750 dilution), [ab219803](#) (1/250 dilution), [ab251613](#) (1/750 dilution), [ab264485](#) (0.5 µg/ml), [ab225894](#) (1/1250 dilution), [ab251596](#) (1/1500 dilution) and [ab96048](#) (10 µg/ml); each using a separate fluorescent tyramide signal amplification system. EDTA based antigen retrieval (Leica Biosystems BOND® Epitope Retrieval Solution 2, pH 9.0, 20 minutes) was used in between rounds of tyramide signal amplification to remove the antibody from

the previous round, to avoid any cross-reactivity. DAPI (dark blue) was used as a nuclear counter stain.

Microscopy and pseudocoloring of individual Opal™ dyes was performed using a Vectra 3 Imaging System (Akoya Biosciences®).

This data is courtesy of ImmunoAtlas and it can be found [here](#).

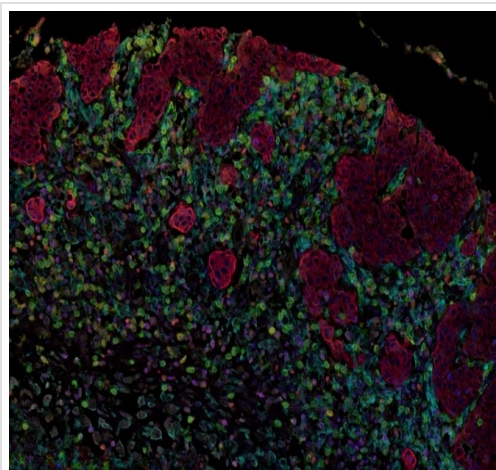

Multiplex immunohistochemistry - Anti-Granzyme B antibody [EPR20129-217] (ab208586)

This image is courtesy of ImmunoAtlas.

Fluorescence multiplex immunohistochemical analysis of Human breast cancer tissue (formalin-fixed paraffin-embedded section).


Merged staining of Anti-PD-L1 ([ab251611](#); cyan; Opal™ 520), Anti-Granzyme B ([ab219803](#); yellow; Opal™ 540), Anti-PD1 ([ab251613](#); magenta; Opal™ 570), Anti-pan Cytokeratin ([ab264485](#); red; Opal™ 620), Anti-EpCAM ([ab225894](#); red; Opal™ 620), Anti-CD8 alpha ([ab251596](#); green; Opal™ 650) and Anti-FOXP3 ([ab96048](#); orange; Opal™ 690). EpCAM and pan-cytokeratin share the same dye and color. Dyes are pseudo-colored for better contrast of the markers.

The immunostaining was performed on a Leica Biosystems BOND® MAX instrument with an Opal™ 6-Plex Detection Kit (NEL821001KT, Akoya Biosciences®).

The section was incubated in six rounds of staining; sequentially for [ab251611](#) (1/750 dilution), [ab219803](#) (1/250 dilution), [ab251613](#) (1/750 dilution), [ab264485](#) (0.5 µg/ml), [ab225894](#) (1/1250 dilution), [ab251596](#) (1/1500 dilution) and [ab96048](#) (10 µg/ml); each using a separate fluorescent tyramide signal amplification system. EDTA based antigen retrieval (Leica Biosystems BOND® Epitope Retrieval Solution 2, pH 9.0, 20 minutes) was used in between rounds of tyramide signal amplification to remove the antibody from the previous round, to avoid any cross-reactivity. DAPI (dark blue) was used as a nuclear counter stain.

Microscopy and pseudocoloring of individual Opal™ dyes was performed using a Vectra 3 Imaging System (Akoya Biosciences®).

This data is courtesy of ImmunoAtlas and it can be found [here](#).

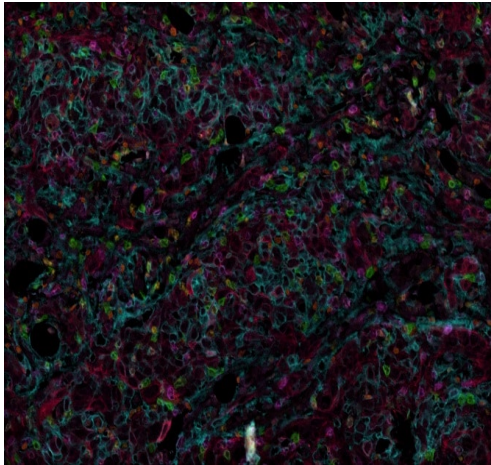

Multiplex immunohistochemistry - Anti-Granzyme B antibody [EPR20129-217] (ab208586)

This image is courtesy of ImmunoAtlas.

Fluorescence multiplex immunohistochemical analysis of Human breast cancer tissue (formalin-fixed paraffin-embedded section).


Merged staining of Anti-PD-L1 ([ab251611](#); cyan; Opal™ 520), Anti-Granzyme B ([ab219803](#); yellow; Opal™ 540), Anti-PD1 ([ab251613](#); magenta; Opal™ 570), Anti-pan Cytokeratin ([ab264485](#); red; Opal™ 620), Anti-EpCAM ([ab225894](#); red; Opal™ 620), Anti-CD8 alpha ([ab251596](#); green; Opal™ 650) and Anti-FOXP3 ([ab96048](#); orange; Opal™ 690). EpCAM and pan-cytokeratin share the same dye and color. Dyes are pseudo-colored for better contrast of the markers.

The immunostaining was performed on a Leica Biosystems BOND® MAX instrument with an Opal™ 6-Plex Detection Kit (NEL821001KT, Akoya Biosciences®).

The section was incubated in six rounds of staining; sequentially for [ab251611](#) (1/750 dilution), [ab219803](#) (1/250 dilution), [ab251613](#) (1/750 dilution), [ab264485](#) (0.5 µg/ml), [ab225894](#) (1/1250 dilution), [ab251596](#) (1/1500 dilution) and [ab96048](#) (10 µg/ml); each using a separate fluorescent tyramide signal amplification system. EDTA based antigen retrieval (Leica Biosystems BOND® Epitope Retrieval Solution 2, pH 9.0, 20 minutes) was used in between rounds of tyramide signal amplification to remove the antibody from the previous round, to avoid any cross-reactivity. DAPI (dark blue) was used as a nuclear counter stain.

Microscopy and pseudocoloring of individual Opal™ dyes was performed using a Vectra 3 Imaging System (Akoya Biosciences®).

This data is courtesy of ImmunoAtlas and it can be found [here](#).

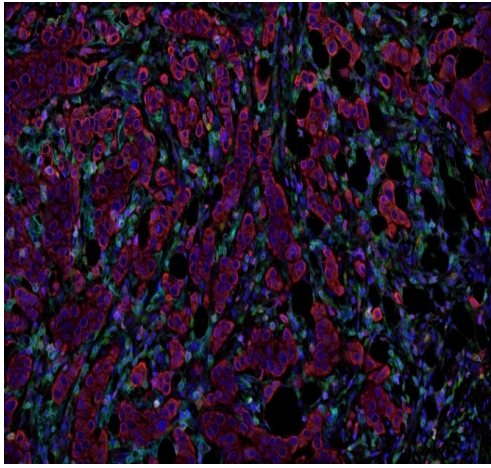

Multiplex immunohistochemistry - Anti-Granzyme B antibody [EPR20129-217] (ab208586)

This image is courtesy of ImmunoAtlas.

Fluorescence multiplex immunohistochemical analysis of Human breast cancer tissue (formalin-fixed paraffin-embedded section).


Merged staining of Anti-PD-L1 ([ab251611](#); cyan; Opal™ 520), Anti-Granzyme B ([ab219803](#); yellow; Opal™ 540), Anti-PD1 ([ab251613](#); magenta; Opal™ 570), Anti-pan Cytokeratin ([ab264485](#); red; Opal™ 620), Anti-EpCAM ([ab225894](#); red; Opal™ 620), Anti-CD8 alpha ([ab251596](#); green; Opal™ 650) and Anti-FOXP3 ([ab96048](#); orange; Opal™ 690). EpCAM and pan-cytokeratin share the same dye and color. Dyes are pseudo-colored for better contrast of the markers.

The immunostaining was performed on a Leica Biosystems BOND® MAX instrument with an Opal™ 6-Plex Detection Kit (NEL821001KT, Akoya Biosciences®).

The section was incubated in six rounds of staining; sequentially for [ab251611](#) (1/750 dilution), [ab219803](#) (1/250 dilution), [ab251613](#) (1/750 dilution), [ab264485](#) (0.5 µg/ml), [ab225894](#) (1/1250 dilution), [ab251596](#) (1/1500 dilution) and [ab96048](#) (10 µg/ml); each using a separate fluorescent tyramide signal amplification system. EDTA based antigen retrieval (Leica Biosystems BOND® Epitope Retrieval Solution 2, pH 9.0, 20 minutes) was used in between rounds of tyramide signal amplification to remove the antibody from the previous round, to avoid any cross-reactivity. DAPI (dark blue) was used as a nuclear counter stain.

Microscopy and pseudocoloring of individual Opal™ dyes was performed using a Vectra 3 Imaging System (Akoya Biosciences®).

This data is courtesy of ImmunoAtlas and it can be found [here](#).

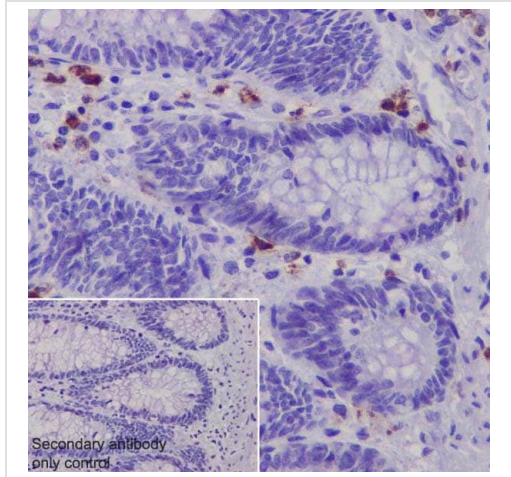

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Granzyme B antibody [EPR20129-217] (ab208586)

Immunohistochemical analysis of paraffin-embedded human colon tissue labeling Granzyme B with ab208586 at 1/250 dilution, followed by Goat Anti-Rabbit IgG H&L (HRP) ([ab97051](#)) at 1/500 dilution. Cytoplasmic staining on some stromal cells of human colon is observed. Counter stained with Hematoxylin.

Secondary antibody only control: Used PBS instead of primary antibody, secondary antibody is Goat Anti-Rabbit IgG H&L (HRP) ([ab97051](#)) at 1/500 dilution.

Perform heat mediated antigen retrieval with Tris/EDTA buffer pH 9.0 before commencing with IHC staining protocol.

Tissue Microarray (TMA) data for ab208586

Normal tissue samples				Malignant tissue samples			
Human cardiac muscle	x	Human placenta	x	Clear cell carcinoma of human kidney	x (immune cells ✓)	Human glioma	x
Human cerebrum	x	Human skeletal muscle	x	Human bladder cancer	x	Human hepatocellular carcinoma	x (immune cells ✓)
Human colon	x (immune cells ✓)	Human skin	x	Human breast carcinoma	x (immune cells ✓)	Human lung carcinoma	x
Human endometrium	x	Human spleen	✓	Human cervical carcinoma	x (immune cells ✓)	Human ovarian carcinoma	x (immune cells ✓)
Human kidney	x	Human stomach	x	Human colon carcinoma	x (immune cells ✓)	Human pancreatic carcinoma	x (immune cells ✓)
Human liver	x	Human testis	x	Human endometrial carcinoma	x (immune cells ✓)	Human prostatic hyperplasia	x (immune cells ✓)
Human lung	x	Human thyroid	x	Human gastric adenocarcinoma	x (immune cells ✓)	Human thyroid carcinoma	x (immune cells ✓)
Human mammary gland	x	Human tonsil	✓				
Human pancreas	x (immune cells ✓)						

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-Granzyme B antibody [EPR20129-217] (ab208586)

Tissue Microarrays stained for "Anti-Granzyme B antibody [EPR20129-217]" using "ab208586" in immunohistochemical analysis. This table provides a detailed overview of positive (tick mark) and negative (cross mark) staining per sample type tested. The sections were pre-treated using Heat mediated antigen retrieval using Bond™ Epitope Retrieval Solution 2 (pH 9.0) for 20 minutes. The sections were incubated with ab208586 for 30 mins at room temperature followed by a ready to use Rabbit specific IHC polymer detection kit HRP/DAB ([ab209101](#)). The immunostaining was performed on a Leica Biosystems BOND® RX instrument.

Why choose a recombinant antibody?



Research with confidence
Consistent and reproducible results



Long-term and scalable supply
Recombinant technology



Success from the first experiment
Confirmed specificity



Ethical standards compliant
Animal-free production

Anti-Granzyme B antibody [EPR20129-217]
(ab208586)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.cn/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors