abcam

Product datasheet

Anti-PD-L1 antibody [SP142] - C-terminal, prediluted ab228463



重组 RabMAb

21 References 11 图像

概述

产品名称 Anti-PD-L1抗体[SP142] - C-terminal, prediluted

描述 兔单克隆抗体[SP142] to PD-L1 - C-terminal, prediluted

宿主 Rabbit

适用于: IHC-P 经测试应用 种属反应性 与反应: Human

免疫原 Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.

阳性对照 IHC-P: Human placenta, tonsil, lung squamous cell carcinoma, cervical squamous cell carcinoma,

skin squamous cell carcinoma, Hodgkin's lymphoma, pancreatic adenocarcinoma and prostate

adenocarcinoma tissues.

常规说明 This product is FOR RESEARCH USE ONLY. For commercial use, please contact

partnerships@abcam.com.

性能

形式 Liquid

存放说明 Shipped at 4°C. Store at +4°C. Do Not Freeze.

存储溶液

Preservative: 0.1% Sodium azide

Constituents: Tris buffered saline, 1% BSA

纯度 Protein A/G purified

纯化说明 Purified from TCS by protein A/G.

克隆 单克隆 克隆编号 SP142 同种型 ΙgG

应用

Abpromise™承诺保证使用ab228463于以下的经测试应用 The Abpromise guarantee


"应用说明"部分 下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。

应 用	Ab评论	说明
IHC-P		1/1. Ready to use. Boil tissue section in EDTA buffer, pH 8.0 for 10 minutes followed by cooling at room temperature for 20 minutes. Primary antibody incubation for 10 minutes at room temperature.

靶标


功能	Involved in the costimulatory signal, essential for T-cell proliferation and production of IL10 and IFNG, in an IL2-dependent and a PDCD1-independent manner. Interaction with PDCD1 inhibits T-cell proliferation and cytokine production.
组织 特异性	Highly expressed in the heart, skeletal muscle, placenta and lung. Weakly expressed in the thymus, spleen, kidney and liver. Expressed on activated T- and B-cells, dendritic cells, keratinocytes and monocytes.
序列相似性	Belongs to the immunoglobulin superfamily. BTN/MOG family. Contains 1 lg-like C2-type (immunoglobulin-like) domain. Contains 1 lg-like V-type (immunoglobulin-like) domain.
细胞定位	Cell membrane and Endomembrane system.

图片

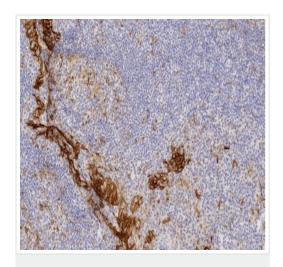

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] - C-terminal, prediluted (ab228463)


IHC image of <u>ab228462</u> staining PD-L1 in human tonsil formalin fixed paraffin embedded tissue sections*, performed on a Leica BOND RX (standard Protocol F, Polymer Refine kit). The section was pre-treated using heat mediated antigen retrieval with EDTA buffer (pH9, epitope retrieval solution 2) for 30 mins at 98°C. The section was then incubated with <u>ab228462</u>, 1/400 working dilution, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system for 8 minutes at room temperature. DAB was used as the chromogen for 10 minutes at room temperature. The section was then counterstained with hematoxylin, blued, dehydrated, cleared and mounted with DPX.

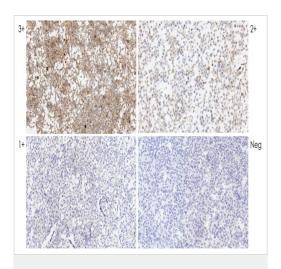

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] - C-terminal, prediluted (ab228463)


IHC image of <u>ab228462</u> staining PD-L1 in PD-L1 Dynamic Range Analyte Control formalin fixed paraffin embedded cell lines (<u>HistoCyte Laboratories</u>), performed on a Leica BOND RX (standard Protocol F, Polymer Refine kit). The section was pretreated using heat mediated antigen retrieval with EDTA buffer (pH9, epitope retrieval solution 2) for 30 mins at 98°C. The section was then incubated with <u>ab228462</u>, 1/400 working dilution, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system for 8 minutes at room temperature. DAB was used as the chromogen for 10 minutes at room temperature. The section was then counterstained with hematoxylin, blued, dehydrated, cleared and mounted with DPX.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

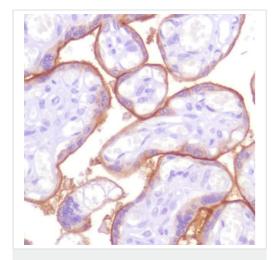


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] - C-terminal, prediluted (ab228463)

Formalin-fixed, paraffin-embedded human placenta tissue stained for PD-L1 using ab228463 in immunohistochemical analysis.

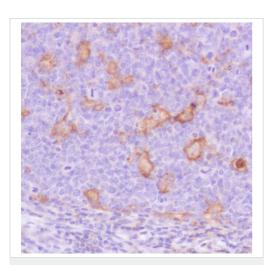


Formalin-fixed, paraffin-embedded human tonsil tissue stained for PD-L1 using ab228463 in immunohistochemical analysis.

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] - C-terminal, prediluted (ab228463)

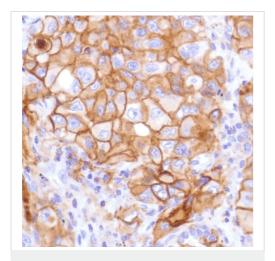


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] - C-terminal, prediluted (ab228463)

Formalin-fixed, paraffin-embedded human lung squamous cell carcinoma tissue stained for PD-L1 using ab228463 in immunohistochemical analysis.

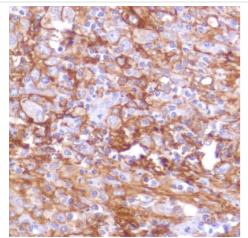


Formalin-fixed, paraffin-embedded human Hodgkin's lymphoma tissue stained for PD-L1 using ab228463 in immunohistochemical analysis.



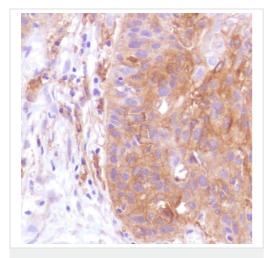


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] -C-terminal, prediluted (ab228463)

Formalin-fixed, paraffin-embedded human cervical squamous cell carcinoma tissue stained for PD-L1 using ab228463 in immunohistochemical analysis.

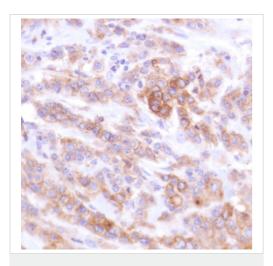


Formalin-fixed, paraffin-embedded human pancreatic adenocarcinoma tissue stained for PD-L1 using ab228463 in immunohistochemical analysis.

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] - C-terminal, prediluted (ab228463)

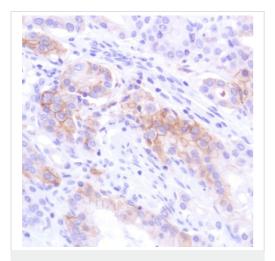


Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] - C-terminal, prediluted (ab228463)

Formalin-fixed, paraffin-embedded human prostate adenocarcinoma tissue stained for PD-L1 using ab228463 in immunohistochemical analysis.

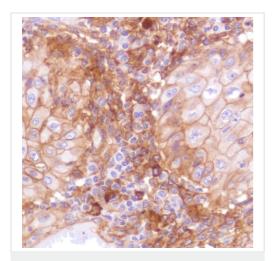

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-PD-L1 antibody [SP142] - C-terminal, prediluted (ab228463)

Formalin-fixed, paraffin-embedded human skin squamous cell carcinoma tissue stained for PD-L1 using ab228463 in immunohistochemical analysis.

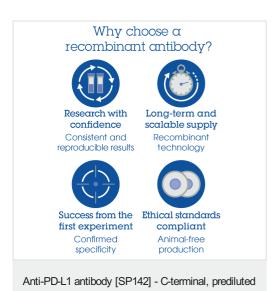

(ab228463)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish

- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.cn/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors