abcam

Product datasheet

Anti-Werner's syndrome helicase WRN antibody ab200

★★★★★ 4 Abreviews 32 References 1 图像

概述

产**品名称** Anti-Werner's syndrome helicase WRN抗体

描述 兔多克隆抗体to Werner's syndrome helicase WRN

宿主 Rabbit

经测试应用 适用于: IHC-P **种属反应性 与反应:** Human

免疫原 Synthetic peptide corresponding to Human Werner's syndrome helicase WRN aa 1-550.

常规说明

The Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

性能

形式 Liquid

存放说明 Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C or -

80°C. Avoid freeze / thaw cycle.

存储溶液 Preservative: 0.05% Sodium azide

纯**度** Whole antiserum

应用

The Abpromise guarantee Abpromise™承诺保证使用ab200于以下的经测试应用

"应用说明"部分下显示的仅为推荐的起始稀释度;实际最佳的稀释度/浓度应由使用者检定。

1

应 用	Ab评论	说明
IHC-P		Use at an assay dependent concentration.

靶标

功能

Multifunctional enzyme that has both magnesium and ATP-dependent DNA-helicase activity and 3'->5' exonuclease activity towards double-stranded DNA with a 5'-overhang. Has no nuclease activity towards single-stranded DNA or blunt-ended double-stranded DNA. Binds preferentially to DNA substrates containing alternate secondary structures, such as replication forks and Holliday junctions. May play an important role in the dissociation of joint DNA molecules that can arise as products of homologous recombination, at stalled replication forks or during DNA repair. Alleviates stalling of DNA polymerases at the site of DNA lesions. Important for genomic integrity. Plays a role in the formation of DNA replication focal centers; stably associates with foci elements generating binding sites for RP-A.

疾病相关

Defects in WRN are a cause of Werner syndrome (WRN) [MIM:277700]. WRN is a rare autosomal recessive progeroid syndrome characterized by the premature onset of multiple agerelated disorders, including atherosclerosis, cancer, non-insulin-dependent diabetes mellitus, ocular cataracts and osteoporosis. The major cause of death, at a median age of 47, is myocardial infarction. Currently all known WS mutations produces prematurely terminated proteins.

Defects in WRN may be a cause of colorectal cancer (CRC) [MIM:114500].

序列相似性

Belongs to the helicase family. RecQ subfamily.

Contains 1 3'-5' exonuclease domain.

Contains 1 helicase ATP-binding domain.

Contains 1 helicase C-terminal domain.

Contains 1 HRDC domain.


翻译后修饰

Phosphorylated by PRKDC. Phosphorylated upon DNA damage, probably by ATM or ATR.

细胞定位

Nucleus > nucleolus. Nucleus.

图片

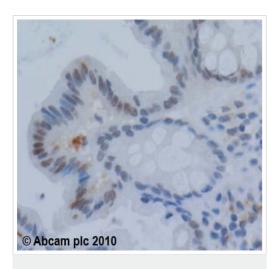

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-Werner's syndrome helicase WRN antibody (ab200)

ab200 (4µg/ml) staining Werner's syndrome helicase WRN in human colon using an automated system (DAKO Autostainer Plus). Using this protocol there is nuclear and some cytoplasmic staining. Sections were rehydrated and antigen retrieved with the Dako 3 in 1 AR buffer citrate pH 6.0 in a DAKO PT link. Slides were peroxidase blocked in 3% H2O2 in methanol for 10 mins. They were then blocked with Dako Protein block for 10 minutes (containing casein 0.25% in PBS) then incubated with primary antibody for 20 min and detected with Dako Envision Flex amplification kit for 30 minutes. Colorimetric detection was completed with Diaminobenzidine for 5 minutes. Slides were counterstained with Haematoxylin and coverslipped under DePeX. Please note that, for manual staining, optimization of primary antibody concentration and incubation time is recommended. Signal amplification may be required.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- · Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.cn/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors